

Maestro2

Robotic 3D Optical Coherence Tomography
with True Color Fundus Camera

**Simple Capture,
Sophisticated Analyses,
Clinical Confidence**



OCTA
AVAILABLE



Maestro2



The Iconic Robotic OCT Fundus Camera, Powered with OCTA

Widefield OCT Scan

12x9mm widefield OCT scan encompasses both macula and disc with thickness metrics and reference database for a comprehensive assessment of eye health.

3D Wide Trend Analysis

Comprehensive longitudinal assessment of optic nerve head photographs, RNFL and ganglion cell layer thickness data in a change-over-time bilateral report.

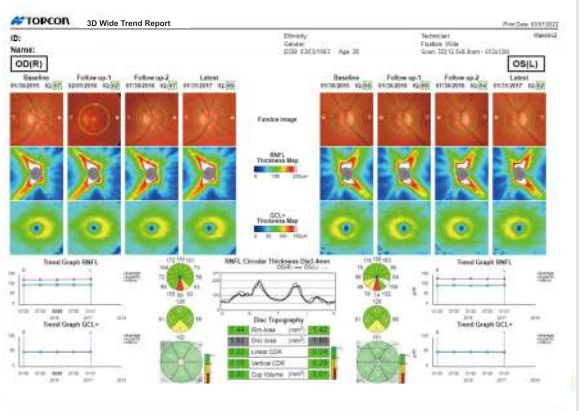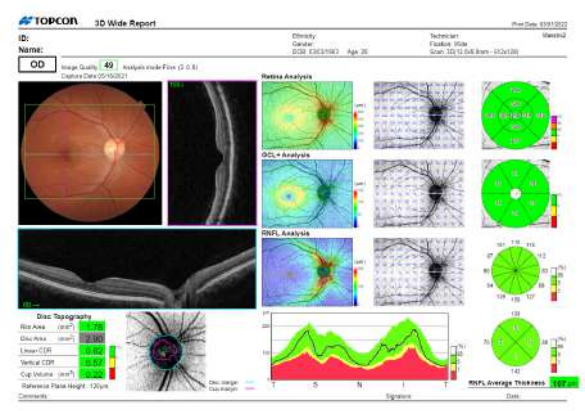

Overview



User-friendly
Robotic OCT +
Fundus Camera



True Color¹
Fundus Photography



Full 360°
Rotating Touchscreen



12x9mm
3D Wide Scan



Anterior Segment OCT



Reference Database



Single Touch Capture




OCTA²

¹ True, full color fundus image simultaneously captured with white light, 24-bit color
² OCTA feature is optional

OCTA


Maestro2 OCTA provides visualization of vascular structure, without the need for contrast dye, using en-face angiographic and cross-sectional imaging of both normal and abnormal vasculature in the retina and choroid.

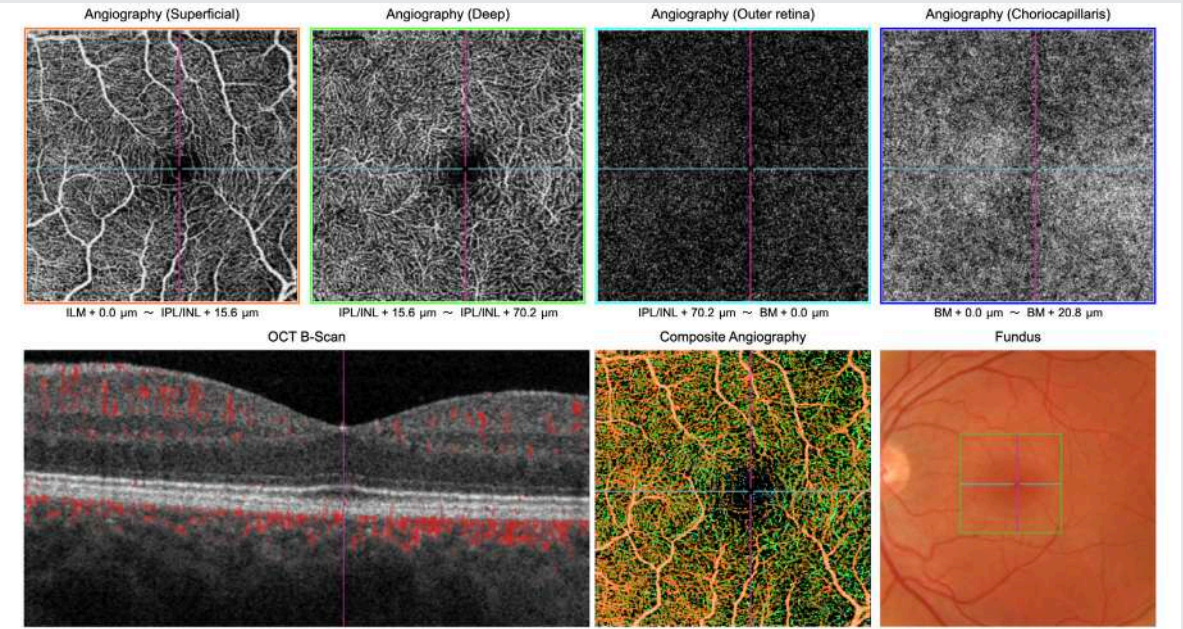


Guidance for Diagnosis

Reports | Retina

Maestro2 provides rich analyses of the macular and optic nerve head regions. Reports can be auto-exported, quickly printed or sent to your image management system or electronic health record(EHR) in common file formats.

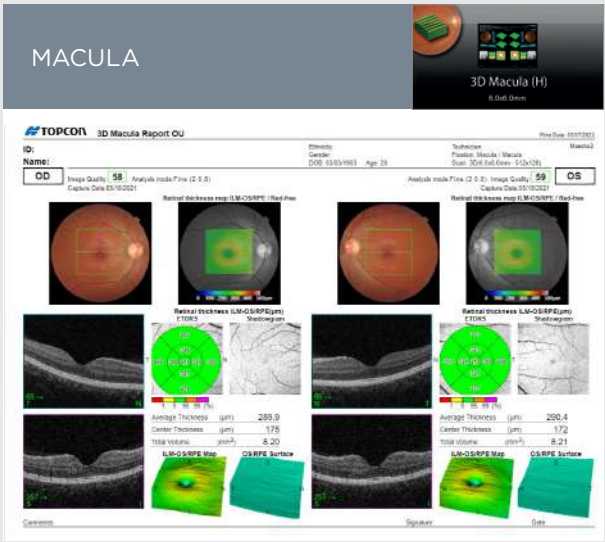


3D MACULA REPORT (OU)
6x6mm bilateral scan report includes true-color and red-free fundus photography with OCT thickness overlay, retinal thickness with reference database, high-resolution OCT scans and thickness surface map

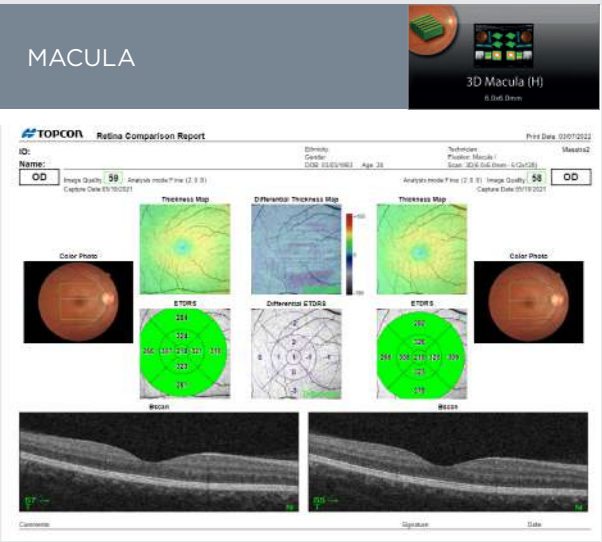

COMPARISON REPORT - CHANGE ANALYSIS
Unilateral visit-to-visit change report with 45° true-color fundus photography, intervisit-registered OCT scans (3D Macula or 3D Wide) and ETDRS. Map and Differential ETDRS displays thickness variance in +/- microns




SINGLE LINE SCAN
45° true-color and red-free fundus photographs with high resolution OCT scan

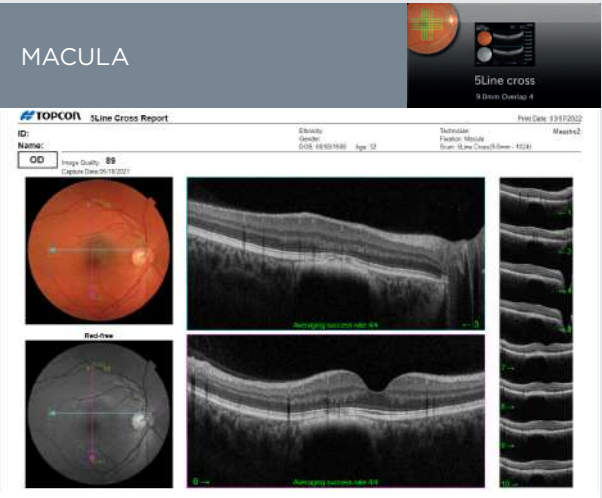


5 LINE CROSS SCAN
5-line cross scan displays horizontal and vertical B-scans (6mm, 9mm)


Reports | Glaucoma

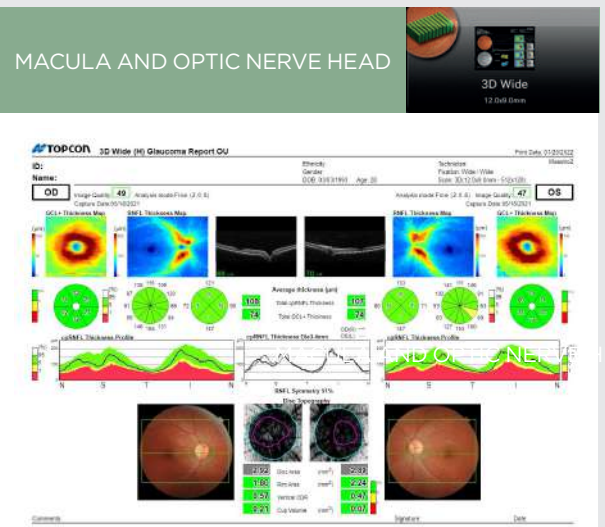


3D WIDE (H) GLAUCOMA REPORT OU
Wide, 12x9mm OU OCT scan report. Includes 45° true-color fundus photography, RNFL thickness, optic nerve head topography, GCL+ thickness; all with reference data

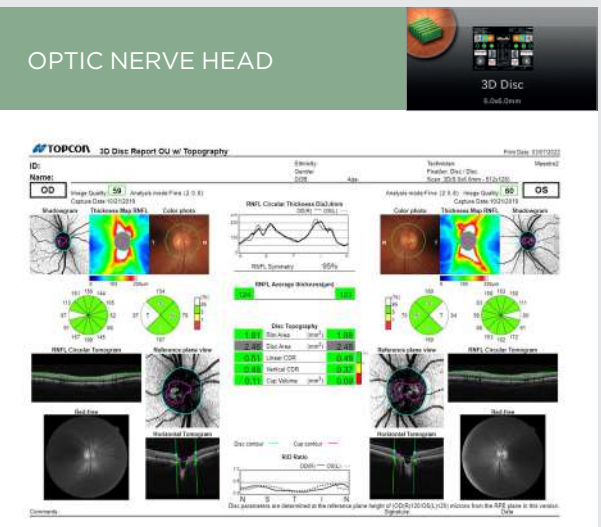

3D DISC REPORT (OU) WITH TOPOGRAPHY
Optic nerve head 6x6mm OCT scans offering conventional analyses with photography in a bilateral report

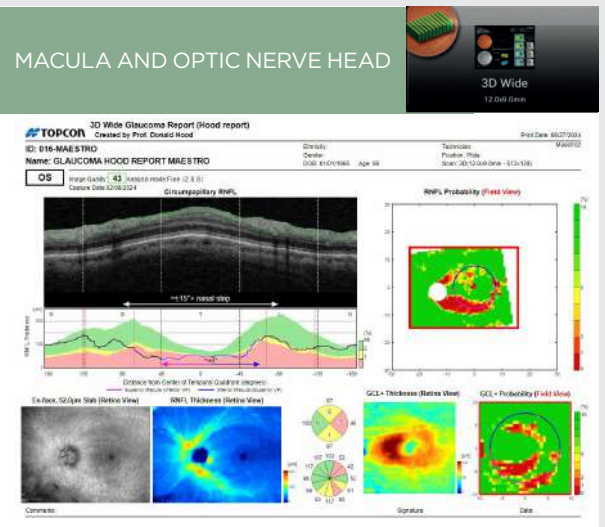

HOOD REPORT FOR GLAUCOMA
Hood Report for Glaucoma featuring probability maps for RNFL and GCL streamlines clinical decision-making

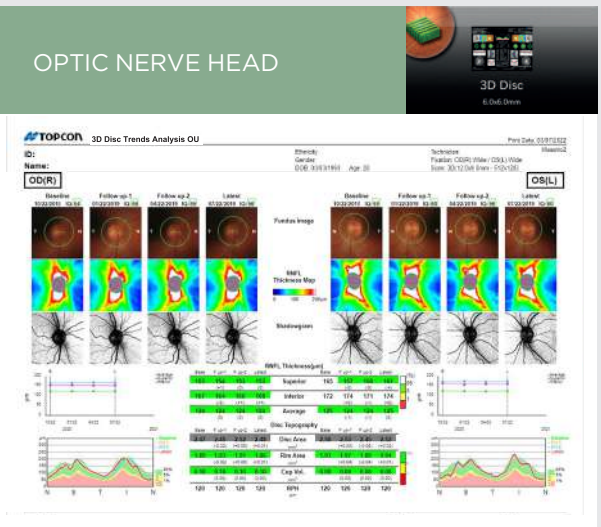

3D DISC TREND ANALYSIS OU
Optic nerve head photographs, RNFL and optic nerve head data presented in a bilateral report show change over time

True Color Fundus Photography¹

Integrated true color fundus camera enables simultaneous capture of the OCT image and fundus photo. PinPoint™ Registration allows multimodal observation of suspected pathology. Small pupil mode and fundus-only capture are also available.

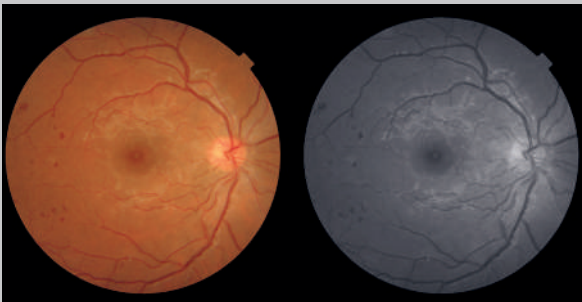


Image courtesy: Michael H. Chen, O.D.

Peripheral Fundus Photography

Automatically select up to nine standard fields or manually manipulate the patient's fixation to create a mosaic image with the IMAGEnet®6 AutoMosaic software.


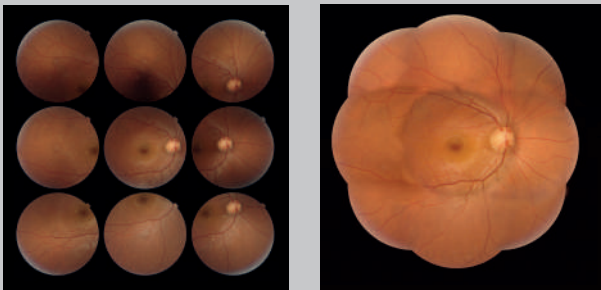


Image courtesy: Michael H. Chen, O.D.


Live Fundus View™ (LFV)

OCT-LFV is a live projection image of the retina that makes the disc, retinal vessels and scanning position easy to see.

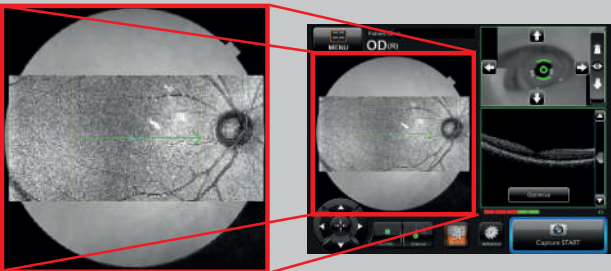


Cataract Mode

Cataract mode automatically adjusts the scanning position to minimize the impact of any opacities such as cataracts.

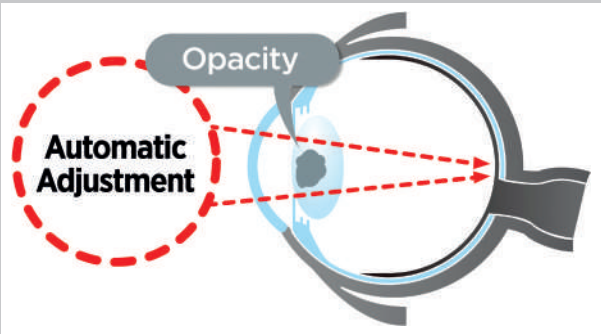

OCTA²

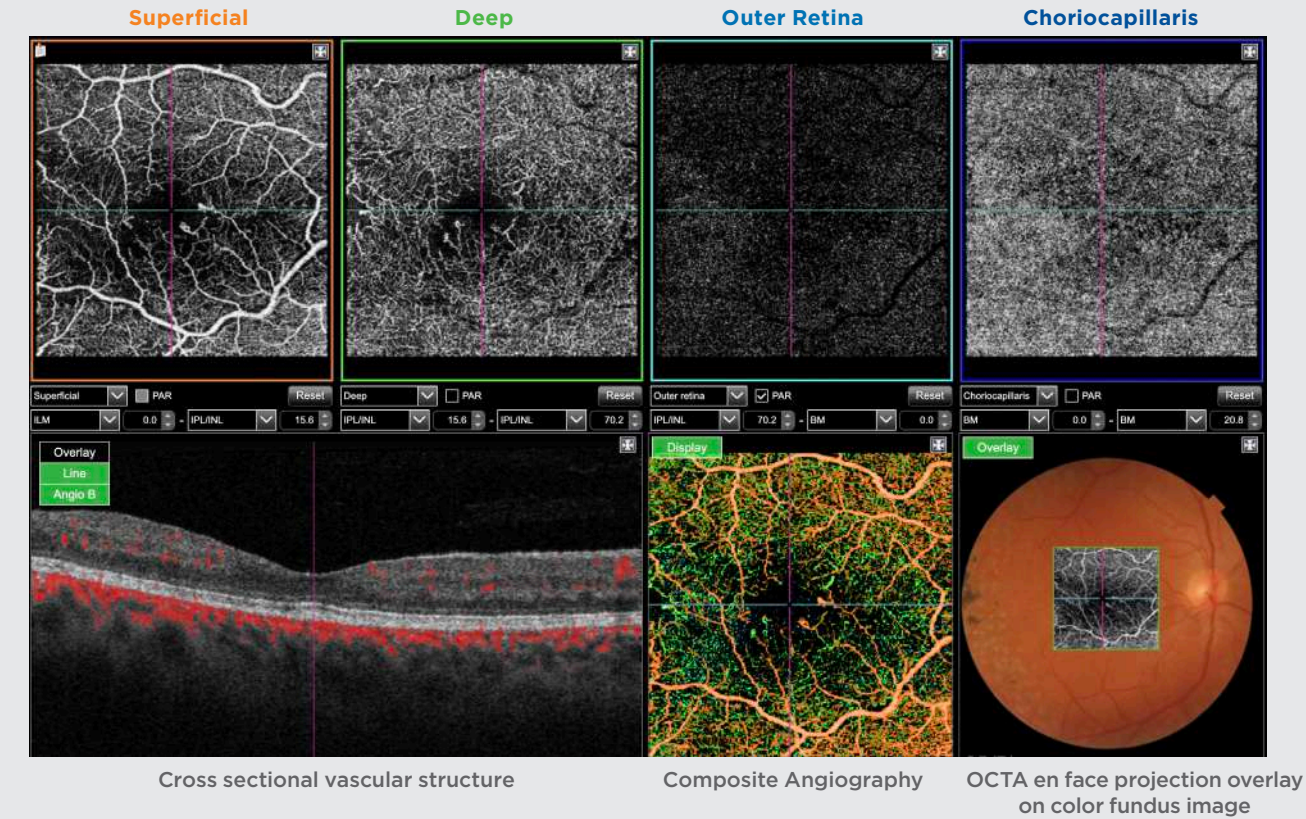


2 OCTA feature is optional

Anterior Segment OCT

Anterior Segment OCT

Capture cornea and anterior chamber angle scan, measure corneal thickness and scleral contact lenses clearance with manual caliper tools.

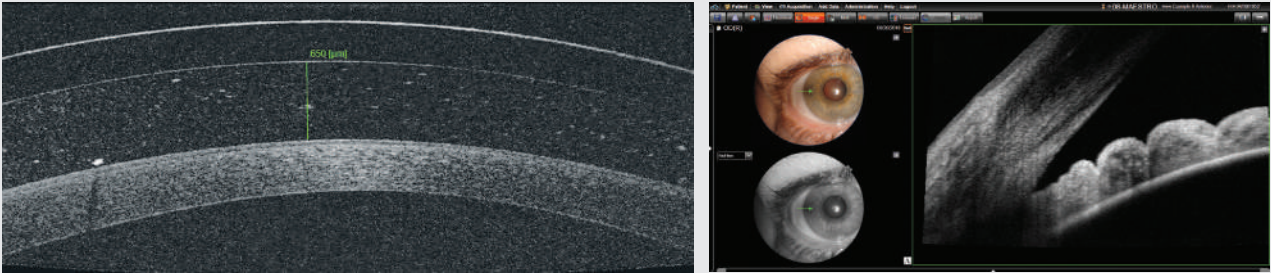

Specifications

Observation & photography of the fundus

Type of photography	Color, Red-free ^(Note 1) & IR ^(Note 3)
Picture angle for photography	45° ± 5% or less 30° or equivalent (digital zoom)
Operating distance	34.8 ± 0.1mm (when taking a picture of fundus)
Photographable diameter of pupil	ø4.0mm or more : When small pupil diaphragm is not used. ø3.3mm or more : When small pupil diaphragm is used.

Observation & photographing of the fundus tomogram

Scan range (on fundus)	Horizontal direction Vertical direction	3 - 12mm ±5% or less 3 - 9mm ±5% or less
OCT scan patterns	3D Wide, 3D Disc, 3D Macula, Line, 5-Line Cross	
OCTA scan size	6 x 6mm, 4.5 x 4.5mm, 3 x 3mm	
Scan speed	50,000 A-Scans per second	
Lateral resolution	20μm or less	
In-depth resolution	6μm or less Pixel spacing: 2.6μm ±2%	
Photographable diameter of pupil	ø2.5mm or more	

Observation & photographing of anterior segment

Type of photography	Color & IR ^(Note 2)
Operating distance	62.6 ± 0.1mm (when taking a picture of anterior segment ^(Note 2))

Observation & photographing of the anterior segment tomogram

Operating distance	62.6 ± 0.1mm (when taking a picture of anterior segment ^(Note 2))	
Scan range (on cornea) ^(Note 3)	Horizontal direction Vertical direction	3 - 6mm ± 5% or less 3 - 6mm ± 5% or less
Scan pattern	Line Anterior Seg, Radial Anterior Seg	
Scan speed	50,000 A-Scans per second	

Electric rating / Dimensions & weight

Source voltage Power input	AC 100 - 240V 50-60Hz 70 - 150VA
Dimensions	13 -19 inches (W) x 21-27 inches (D) x 20-29 inches (H) / 340 - 480 mm (W) x 543 - 680 mm (D) x 500 - 735 mm (H)
Weight	55 lbs / 25kg

(Note 1) Digital red-free photography that processes a color image and displays it in pseudo-red-free condition

(Note 2) When the attachment for anterior segment is included in the system configuration

(Note 3) This is used only for recording the position where a tomogram is captured

* Product name is 3D Optical Coherence Tomography, 3D OCT-1 (Type: Maestro2)

TOPCON HEALTHCARE, INC.

111 Bauer Drive, Oakland, NJ 07436, U.S.A.

Phone: +1-201-599-5100

www.topconhealthcare.com

IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.
Not all products, services, or offers are available in all markets. Contact your local distributor for country-specific information and availability.



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, JAPAN.

Phone: +81-(0)3-3558-2522/2502

Fax: +81-(0)3-3965-6898

www.topconhealthcare.jp

